




PROCESS FLOW DIAGRAM FOR AEROSOL MANUFACTURING

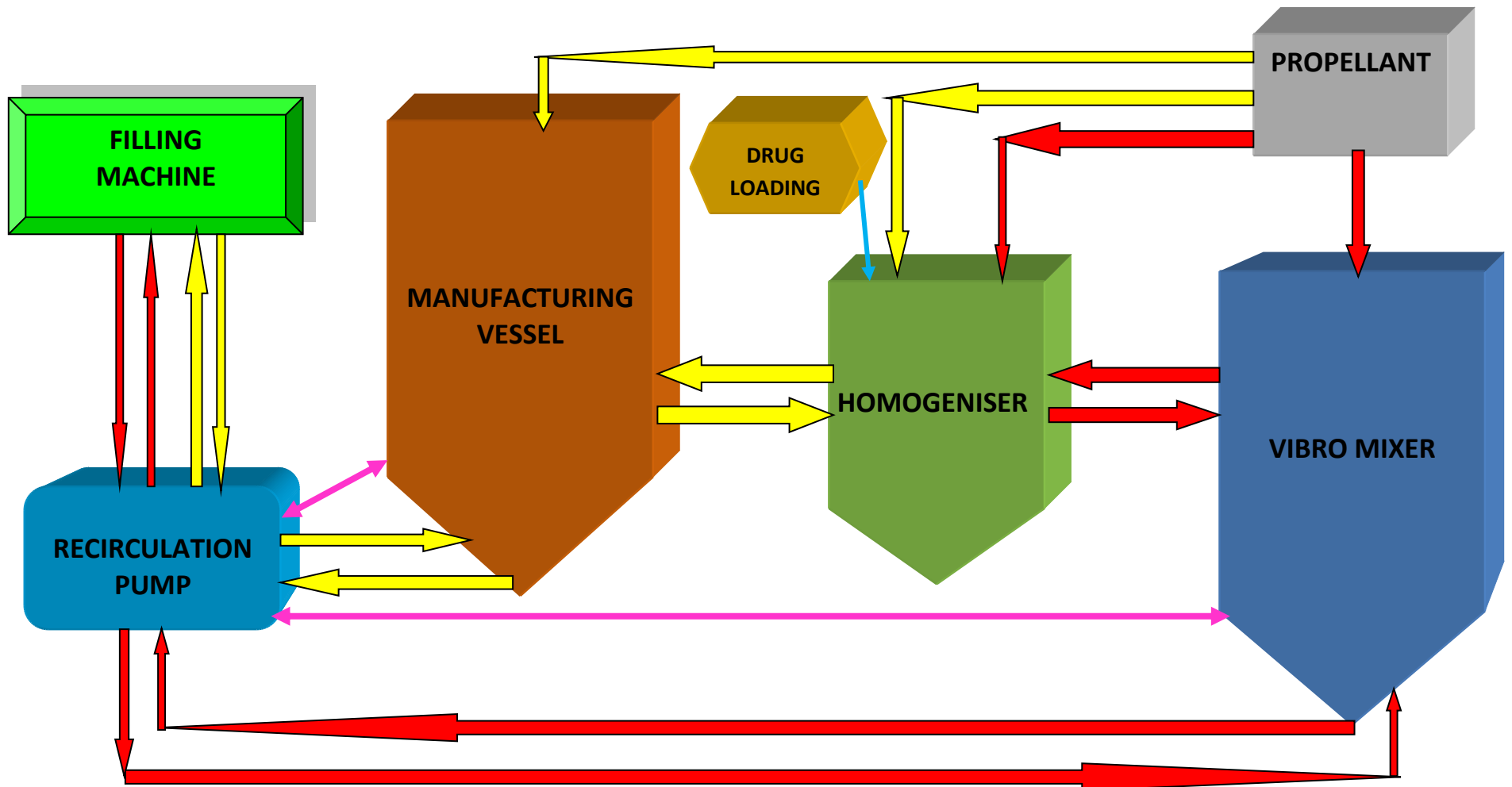

Yellow arrow indicates the process flow when Manufacturing Vessel is used and Vibro mixer is closed.

Red arrow indicates the process flow when Vibro Mixer is used and Manufacturing Vessel is closed.



PHARMA DEVILS

QUALITY ASSURANCE DEPARTMENT

PROCESS FLOW DIAGRAM FOR AEROSOL MANUFACTURING

As per the equipment suitability, and process flow is concerned, both Manufacturing Vessel and Vibro Mixer cannot be operated simultaneously. When Vibro mixer is in use, Propellant inlet, Recirculation Pump, Homogeniser, Return Pipelines from filling machine will be connected with Vibro Mixer and Outlet of Manufacturing vessel will remain Closed and vice versa. Hence there is no chance of product mix-up.